

Thaw and Culture Details

Cell Line Name	JHU172i		
WiCell Lot Number	DB36377		
Provider	Johns Hopkins University – Laboratory of Dr. Lewis Becker		
Banked By	Johns Hopkins University – Laboratory of Dr. Lewis Becker		
Thaw and Culture Recommendations	WiCell recommends thawing 1 vial into 4 wells of a 6 well plate. WiCell recommends thawing using ROCK Inhibitor for best results.		
Culture Platform	Feeder Independent		
	Medium: E8		
	Matrix: Vitronectin		
Protocol	WiCell Feeder Independent E8 Medium Protocol		
Passage Number	p10 These cells were cultured for 10 passages post reprogramming prior to freeze. Add +1 to the passage number to best represent the overall passage number of the cells at thaw.		
Date Vialed	26-January-2016		
Vial Label	P172 P10 1/26/16 1M cells		
Biosafety and Use Information	This cell line is of human origin. Appropriate biosafety precautions should be followed when working with these cells. The end user is responsible for ensuring that the cells are handled and stored in an appropriate manner. WiCell is not responsible for damages or injuries that may result from the use of these cells. Cells distributed by WiCell are intended for research purposes only and are not intended for use in humans.		

Testing Performed by WiCell

Test Description	Test Provider	Test Method	Test Specification	Result
Karyotype by G-banding	WiCell	SOP-CH-003	Expected karyotype	See Report
Post-Thaw Viable Cell Recovery	WiCell	SOP-CH-305	Recoverable attachment after passage	Pass
Identity by STR	UW Translational Research Initiatives in Pathology Laboratory	PowerPlex 16 HS System by Promega	Defines profile	Pass
Sterility	Steris	ST/07	Negative	Pass
Mycoplasma	WiCell	SOP-CH-044	Negative	Pass

Testing Reported by Provider

The Provider stated that some or all of the additional analyses listed below may have been performed for this cell line. For more information, publication and dbGaP links, where available, are provided on the cell line specific web page on the WiCell website.

- Embryoid bodies
- Infinium® Expanded Multi-Ethnic Genotyping Array (MEGAEX)



Approval Date	Quality Assurance Approval
14-July-2016	R/15,2019 X JKG KG Quality Ausurance Signed by Gay, Jenna



Chromosome Analysis Report: 077613


Date Reported: Tuesday, July 23, 2019

Cell Line: JHU172i-DB36377 14899

Passage#: 11

Date of Sample: 7/18/2019 Specimen: Human IPS

Results: 46,XY

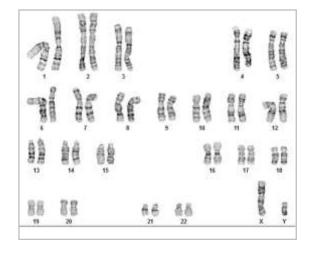

Cell Line Sex: Male

Reason for Testing: lot release testing

Investigator: , WiCell

Cell: 42

Slide: G03

Slide Type: Karyotype

Total Counted: 20 Total Analyzed: 8

Total Karyogrammed: 4

Band Resolution: 475 - 525

Interpretation:

Data.

This is a normal karyotype; no clonal abnormalities were detected at the stated band level of resolution.

Completed by: , CG(ASCP) Reviewed and Interpreted by: , PhD, FACMG

Date:	Sent By:	Sent To:	QC Review By:

Limitations: This assay allows for microscopic visualization of numerical and structural chromosome abnormalities. The size of structural abnormality that can be detected is >3-10Mb, dependent upon the G-band resolution obtained from this specimen. For the purposes of this report, band level is defined as the number of G-bands per haploid genome. It is documented here as "band level", i.e., the range of bands determined from the four karyograms in this assay. Detection of heterogeneity of clonal cell populations in this specimen (i.e.,mosaicism) is limited by the number of metaphase cells examined, documented here as "# of cells counted".

This assay was conducted solely for listed investigator/institution. The results of this assay are for research use only. Unless otherwise mutually agreed in writing, the services provided to you hereunder by WiCell Research Institute, Inc. ("WiCell") are governed solely by WiCell's Terms and Conditions of Service, found at www.wicell.org/privacyandterms. Any terms you may attach to a purchase order or other document that are inconsistent, add to, or conflict with WiCell's Terms and Conditions of Service are null and void and of no legal force or effect.



Short Tandem Repeat HISTOLOGY - IHC - MOLECULAR - IMAGING

WiCell Research Institute

Quality Assurance Department

Requestor:

Analysis



characterization@wicell.org

(608) 316-4145

Receive Date: 07/29/19 **Report Sent:** 08/04/19

Assav Date: 07/30/19

File Name: STR 190731 wmr

Report Date: 08/01/19

TRIP Laboratory (Molecular) https://research.pathology.wisc.edu/trip-home/ (608) 265-9168

Department of Pathology and Laboratory Medicine

Sample Report: 14899-STR

Sample Name on Tube: 14899-STR

93.5 ng/ μ L, (A260/280=2.17)

Sample Type: Cells

Cell Count: ~2 million cells

STR Locus	STR Genotype Repeat #	STR Genotype
FGA	16–18,18.2,19,19.2,20,20.2,21,21.2,22, 22.2, 23, 23.2, 24, 24.2, 25, 25.2, 26–30, 31.2, 43.2, 44.2,45.2, 46.2	Identifying
TPOX	6-13	information has been redacted to
D8S1179	7-18	protect donor
vWA	10-22	confidentiality. If
Amelogenin	X,Y	more information
Penta_D	2.2, 3.2, 5, 7-17	is required,
CSF1PO	6-15	please, contact
D16S539	5, 8-15	WiCell's Technical
D7S820	6-14	Support.
D13S317	7-15	
D5S818	7-16	
Penta_E	5-24	
D18S51	8-10, 10.2, 11-13, 13.2, 14-27	
D21S11	24,24.2,25,25.2,26-28,28.2,29,29.2, 30, 30.2,31, 31.2,32,32.2,33,33.2, 34,34.2,35,35.2,36-38	
TH01	4-9,9.3,10-11,13.3	
D3S1358	12-20	

Results: Based on the 14899-STR cells submitted by WiCell QA dated and received on 07/29/19, this sample (Label on Tube: 14899-STR) defines the STR profile of the human cell line JHU172i comprising 27 allelic polymorphisms across the 15 STR loci analyzed.

Interpretation: No STR polymorphisms other than those corresponding to the human JHU172i cell line were detected and the concentration of DNA required to achieve an acceptable STR genotype (signal/ noise) was equivalent to that required for the standard procedure (~1 ng/amplification reaction) from human genomic DNA. This result suggests that the 14899-STR sample submitted corresponds to the JHU172i cell line and was not contaminated with any other human cells or a significant amount of mouse feeder layer cells.

Sensitivity: Sensitivity limits for detection of STR polymorphisms unique to either this or other human cell lines is $\sim 2-5\%$.

X WMR \mathbf{X} RMB Digitally Signed on 08/04/19 Digitally Signed on 08/04/19 ■ PhD, Director / Co-Director TRIP Laboratory, Molecular UWHC Molecular Diagnostics Laboratory / UWSMPH TRIP Laboratory

Native Product Sterility Report



SAMPLE #:

19071395

WiCell

DATE RECEIVED:

18-Jul-19

504 S Rosa Road, Rm 101

TEST INITIATED:

24-Jul-19

Madison, WI 53719

TEST COMPLETED:

07-Aug-19

SAMPLE NAME / DESCRIPTION:

STAN037i-118-1

DB30906

14904

JHU058i JHU172i DB41092 14905 DB36377

14906

JHU170i JHU225i

DB36371 DB41417

DB36222

14907 14908

JHU143i JHU104i

JHU080i

DB41347 14909 DB41282

14910 14911

14912

JHU097i DB41267 STAN343i-998C1

DB35654

14913

UNIQUE IDENTIFIER:

TEST RESULTS:

	# Positives	
# Tested	(Growth)	- Control
10	0	2 Negatives

TEST SUMMARY:

# Samples	Media Type	Volume (mL)	Incubation Temperature (° C)	Incubation Duration (Days)
10	TSB	40	20-25	14
10	FTG	40	30-35	14

REFERENCE:

Processed according to LAB-003: Sterility Test Procedure

PD #:

000053

TEST METHODOLOGY:

USP - Direct Transfer


COMMENTS:

NA

REVIEWED BY

DATE <u>07 Aug 19</u>

Specific test results may not be indicative of the characteristics of any other samples from the same lot or similar lots. This test report shall not be reproduced, except in full, without prior written approval. Liability is limited to the costs of the tests.Results applied to samples as received.

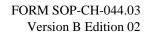

WiCell

Mycoplasma Assay Report

PCR-based assay performed by WiCell WiCell Lot Release Testing 22Jul19

#	Sample Name	Result	Comments/Suggestions
1	JHU172i-DB36377 14899	Negative	Band was not seen at 270bp, indicating the absence of mycoplasma
2	Positive (+) Control	Positive	
3	Negative (-) Control	Negative	

Reported by: Anna Lisa Larson, Senior Research Specialist Reviewed by: Katie Remondini, Cell Culture Specialist

Unless otherwise mutually agreed in writing, the services provided to you hereunder by WiCell Research Institute, Inc. ("WiCell") are governed solely by WiCell's Terms and Conditions of Service, found at www.wicell.org/privacyandterms. Any terms you may attach to a purchase order or other document that are inconsistent, add to, or conflict with WiCell's Terms and Conditions of Service are null and void and of no legal force or effect.

A gel image is available upon request.