



Thaw and Culture Details

Cell Line Name	CREM024i-SS36-1
WiCell Lot Number	WB67440
Provider	Boston University – Laboratory of Dr. Martin Steinberg
Banked By	WiCell
Thaw and Culture Recommendations	WiCell recommends thawing 1 vial into 3 wells of a 6 well plate.
Culture Platform	Feeder Independent
	Medium: mTeSR™Plus
	Matrix: Matrigel®
Protocol	WiCell Feeder Independent mTeSR™Plus Protocol
Passage Number	p12 These cells were cultured for 11 passages prior to freeze and post colony selection. WiCell adds +1 to the passage number at freeze to best represent the overall passage number of the cells at thaw. Plated cells at thaw should be labeled passage 12.
Date Vialied	29-February-2020
Vial Label	CREM024i-SS36-1 p12 WB67440
Biosafety and Use Information	Appropriate biosafety precautions should be followed when working with these cells. The end user is responsible for ensuring that the cells are handled and stored in an appropriate manner. WiCell is not responsible for damages or injuries that may result from the use of these cells. Cells distributed by WiCell are intended for research purposes only and are not intended for use in humans.

Testing Performed by WiCell

Test Description	Test Provider	Test Method	Test Specification	Result
Karyotype by G-banding	WiCell	SOP-CH-003	Expected karyotype	See Report
Post-Thaw Viable Cell Recovery	WiCell	SOP-CH-305	≥ 15 Undifferentiated Colonies prior to passage, ≤ 30% Differentiation prior to passage, and recoverable attachment after passage	Pass
Identity by STR	UW Translational Research Initiatives in Pathology Laboratory	PowerPlex 16 HS System by Promega	Defines STR profile of deposited cell line	Pass
Sterility	Steris	ST/07	Negative	Pass
Mycoplasma	WiCell	SOP-CH-044	Negative	Pass



Testing Reported by Provider

The Provider stated that some or all of the additional analyses listed below may have been performed for this cell line. For more information, publication and dbGaP links, where available, are provided on the cell line specific web page on the WiCell website.

- Digital Genome Sequencing
- Infinium® Expanded Multi-Ethnic Genotyping Array (MEGA^{EX})

Approval Date	Quality Assurance Approval
18-June-2020	<p data-bbox="1019 590 1057 604">6/18/2020</p> <p data-bbox="878 617 927 642">X JKG</p> <p data-bbox="878 646 964 680">JKG Quality Assurance Signed by: Gay, Jenna</p>

Date Reported: Monday, March 23, 2020

Cell Line: CREM024i-SS36-1-WB67440

Passage#: 12

Date of Sample: 3/16/2020


Specimen: Human iPSC

Results: 46,XY

Cell Line Sex: Male

Reason for Testing: LOT_RELEASE

Investigator: WiCell Stem Cell Bank, WiCell

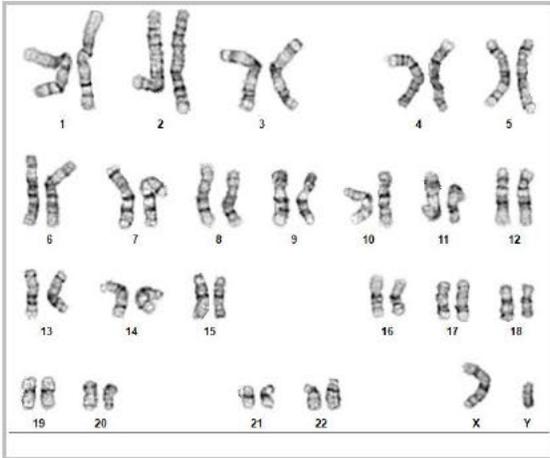

Cell: 41

Slide: G03

Slide Type: Karyotype

Total Counted: 20

Total Analyzed: 8

Total Karyogrammed: 4

Band Resolution: 400 - 450

Interpretation:

This is a normal karyotype; no clonal abnormalities were detected at the stated band level of resolution.

Completed by: [REDACTED], CG(ASCP)

Reviewed and Interpreted by: [REDACTED], Ph.D.

Date: _____ **Sent By:** _____ **Sent To:** _____ **QC Review By:** _____

Limitations: This assay allows for microscopic visualization of numerical and structural chromosome abnormalities. The size of structural abnormality that can be detected is >3-10Mb, dependent upon the G-band resolution obtained from this specimen. For the purposes of this report, band level is defined as the number of G-bands per haploid genome. It is documented here as "band level", i.e., the range of bands determined from the four karyograms in this assay. Detection of heterogeneity of clonal cell populations in this specimen (i.e., mosaicism) is limited by the number of metaphase cells examined, documented here as "# of cells counted".

This assay was conducted solely for listed investigator/institution. The results of this assay are for research use only. Unless otherwise mutually agreed in writing, the services provided to you hereunder by WiCell Research Institute, Inc. ("WiCell") are governed solely by WiCell's Terms and Conditions of Service, found at www.wicell.org/privacyandterms. Any terms you may attach to a purchase order or other document that are inconsistent, add to, or conflict with WiCell's Terms and Conditions of Service are null and void and of no legal force or effect.

Short Tandem Repeat Analysis

Requestor: WiCell Characterization

Receive Date: 03/23/20
 Report Sent: 04/21/20

Label on tube	MCW087i-U7112-WB67434 p.21 (80872)	CREM049i-BR21-1-DB66767 p.16 (80873)	CREM050i-BR23-1-DB66768 p.15 (80874)	WISCe011-A-40-WB67443 p.9 (80875)	SCRPO203i-DB42677 p.11 (80886)	CREM058i-BR43-1-DB66777 p.10 (80895)	CREM054i-BR33-1-DB66773 p.7 (80898)
Label on Report	Identifying information has been redacted to protect donor confidentiality. If more information is required, please, contact WiCell's Technical Support .						
conc (ng/μL)							
A260/280							
Assay Date							
File Name							
FGA							
TPOX							
D8S1179							
vWA							
Amelogenin							
Penta_D							
CSF1PO							
D16S539							
D7S820							
D13S317							
D5S818							
Penta_E							
D18S51							
D21S11							
TH01							
D3S1358							
Allelic Polymorphisms	28	30	26	28	26	28	29
Matches*				80711			
Comments			microvariant				



HISTOLOGY - IHC - MOLECULAR – IMAGING
 Department of Pathology and Laboratory Medicine
 TRIP Laboratory (Molecular)
<https://research.pathology.wisc.edu/trip-home/>
 (608) 265-9168



Your Lab Partner
characterization@wicell.org
 (608) 316-4145

Short Tandem Repeat Analysis

Label on tube	Elf1-WB67433 p.16 (80899)	CREM024i-SS36-1- WB67440 p.12 (80952)	SCRPO302i-DB42682 p.14 (80953)	STAN312i-906C3-DB44421 p.16 (81039)
Label on Report	Identifying information has been redacted to protect donor confidentiality. If more information is required, please, contact WiCell's Technical Support .			
conc (ng/μL)				
A260/280				
Assay Date				
File Name				
FGA				
TPOX				
D8S1179				
vWA				
Amelogenin				
Penta_D				
CSF1PO				
D16S539				
D7S820				
D13S317				
D5S818				
Penta_E				
D18S51				
D21S11				
TH01				
D3S1358				
Allelic Polymorphisms	27	27	26	26
Matches*	16479			
Comments				



HISTOLOGY - IHC - MOLECULAR – IMAGING
 Department of Pathology and Laboratory Medicine
 TRIP Laboratory (Molecular)
<https://research.pathology.wisc.edu/trip-home/>
 (608) 265-9168



Your Lab Partner
 characterization@wicell.org
 (608) 316-4145

Short Tandem Repeat Analysis

Results: Based on the DNA submitted by WiCell Characterization Department dated and received on 03/23/20, these samples define the STR profiles of the human cell lines as indicated by name. The genotypic profiles comprise a range of 26-30 allelic polymorphisms across the 15 STR loci analyzed.

Interpretation: The concentration of DNA required to achieve an acceptable STR genotype (signal/ noise) was equivalent to that required for the standard procedure (~1 ng/amplification reaction) from human genomic DNA. These results suggests that the cells submitted correspond to the cell lines as named and were not contaminated with any other human cells or a significant amount of mouse feeder layer cells.

Sensitivity: Sensitivity limits for detection of STR polymorphisms unique to either this or other human cell lines is ~2-5%.

¹ For sample 80874 a microvariant exists at the D3S1358 loci with a size less than 11 but undefined due to the lack of sizing standard prior to 11 at this loci.

*Acknowledge TRIP in your publications, posters & presentations. For details, see:
<https://research.pathology.wisc.edu/acknowledging-trip/>*

* **Note:** The STR profile of the following sample is an exact match for the given sample/samples.

X *RMB*

Digitally Signed on 04/21/20

██████████, BA

TRIP Laboratory, Molecular

X *WMR*

Digitally Signed on 04/21/20

██████████, PhD, Director / Co-Director

UWHC Molecular Diagnostics Laboratory / UWSMPH TRIP Laboratory

Testing was accomplished by analysis of human genetic polymorphisms at STR loci. This methodology has not yet been approved by the FDA and is for investigational use only.

Unless otherwise mutually agreed in writing, the services provided to you hereunder by WiCell Research Institute, Inc. ("WiCell") are governed solely by WiCell's Terms and Conditions of Service, found at <https://www.wicell.org/media.acux/ca76d97c-862a-43f3-b02a-ab2d1e619100>. Any terms you may attach to a purchase order or other document that are inconsistent, add to, or conflict with WiCell's Terms and Conditions of Service are null and void and of no legal force or effect.



Mycoplasma Assay Report

PCR-based assay performed by WiCell

WiCell

20Mar20

FORM SOP-CH-048.01

Version B Edition 01

Sample Name	Result	Comments/Suggestions
WISCe011-A-40-WB67443 (80870)	Negative	Band was not seen at 270bp, indicating the absence of mycoplasma.
CREM024i-SS36-1-WB67440 (80954)	Negative	Band was not seen at 270bp, indicating the absence of mycoplasma.
CREM048i-BR3-1-DB66766 (81033)	Negative	Band was not seen at 270bp, indicating the absence of mycoplasma.
CREM049i-BR21-1-DB66767 (81034)	Negative	Band was not seen at 270bp, indicating the absence of mycoplasma.
CREM050i-BR23-1-DB66768 (81035)	Negative	Band was not seen at 270bp, indicating the absence of mycoplasma.
CREM054i-BR33-1-DB66773 (81036)	Negative	Band was not seen at 270bp, indicating the absence of mycoplasma.
CREM055i-BR37-1-DB66774 (81037)	Negative	Band was not seen at 270bp, indicating the absence of mycoplasma.
GFAP-R88C-DB67436 (81038)	Negative	Band was not seen at 270bp, indicating the absence of mycoplasma.
Positive (+) Control	Positive	
Negative (-) Control	Negative	

Reported by: [REDACTED], Assistant Cell Culture Specialist

Reviewed by: [REDACTED], Assistant Research Specialist

Unless otherwise mutually agreed in writing, the services provided to you hereunder by WiCell Research Institute, Inc. ("WiCell") are governed solely by WiCell's Terms and Conditions of Service, found at www.wicell.org/privacyandterms. Any terms you may attach to a purchase order or other document that are inconsistent, add to, or conflict with WiCell's Terms and Conditions of Service are null and void and of no legal force or effect.

A gel image is available upon request.